# Gina

Bruksanvisning | User Guide

**tryggitel as**

Tryggitel AS  
Hans Nielsen Hauges gate 50  
0481 Oslo  
Telefon: 22 15 75 90  
firmapost@tryggitel.no  
www.tryggitel.no

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
## Package contents

- Gina carephone
- Erik radio trigger with accessories.
- Telephone lead/telephone adapter.
- User guide

## Carephone parts

There are three buttons on the front, with the following standard functions:

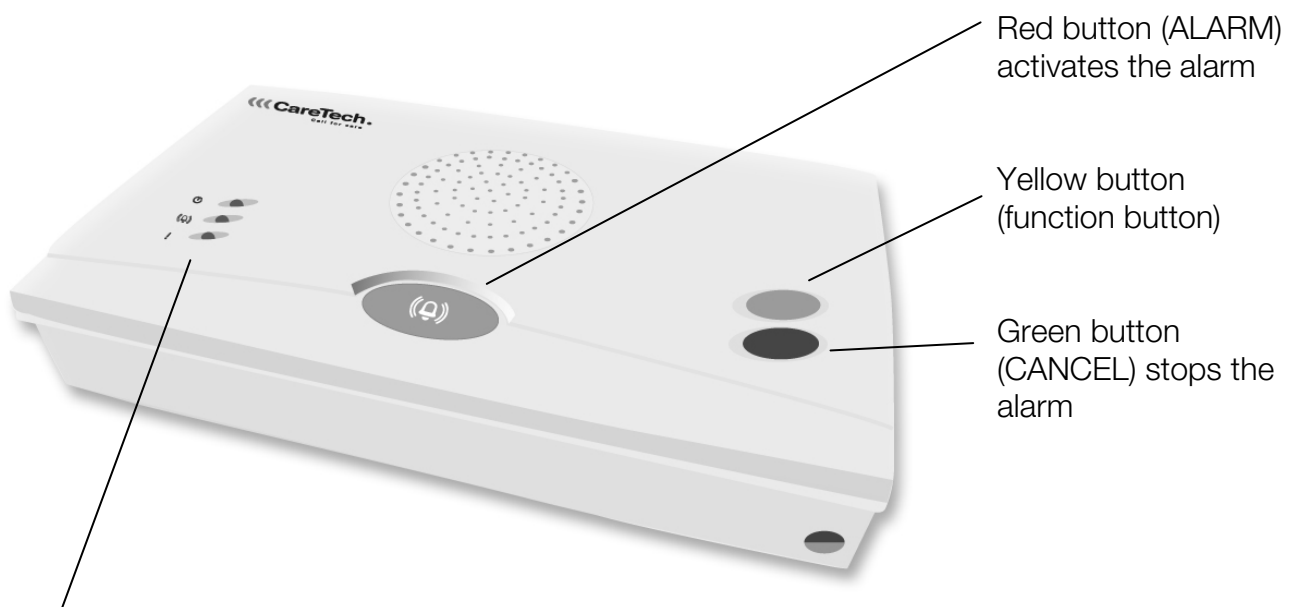The following indicator lights can be found on the front of the carephone:

### **On (Green)** – normally lit during operation

- Steady light: Indicates receipt of radio signals on programming radio alarm.  
Flashing light: Radio programming state.


### **Alarm (Red)**

- Slow flashing: Alarm.  
Rapid flashing: Connection to alarm receiver on alarm.

### **Error (Red)**


- Rapid flashing: No connection to telephone line.  
Slow flashing: Battery fault.

The following connectors can be found on the reverse of the carephone.
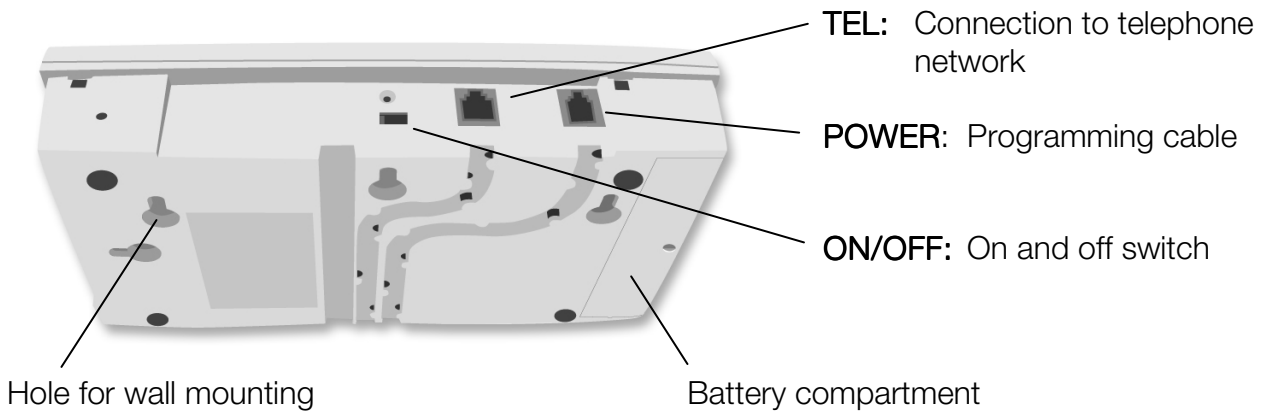
### Radio trigger components

Accessories included: wristband, wristband attachment, neckband, pocket clip, "easy-press" accessory and cover. The accessories are attached to the alarm button with a simple click-in action.
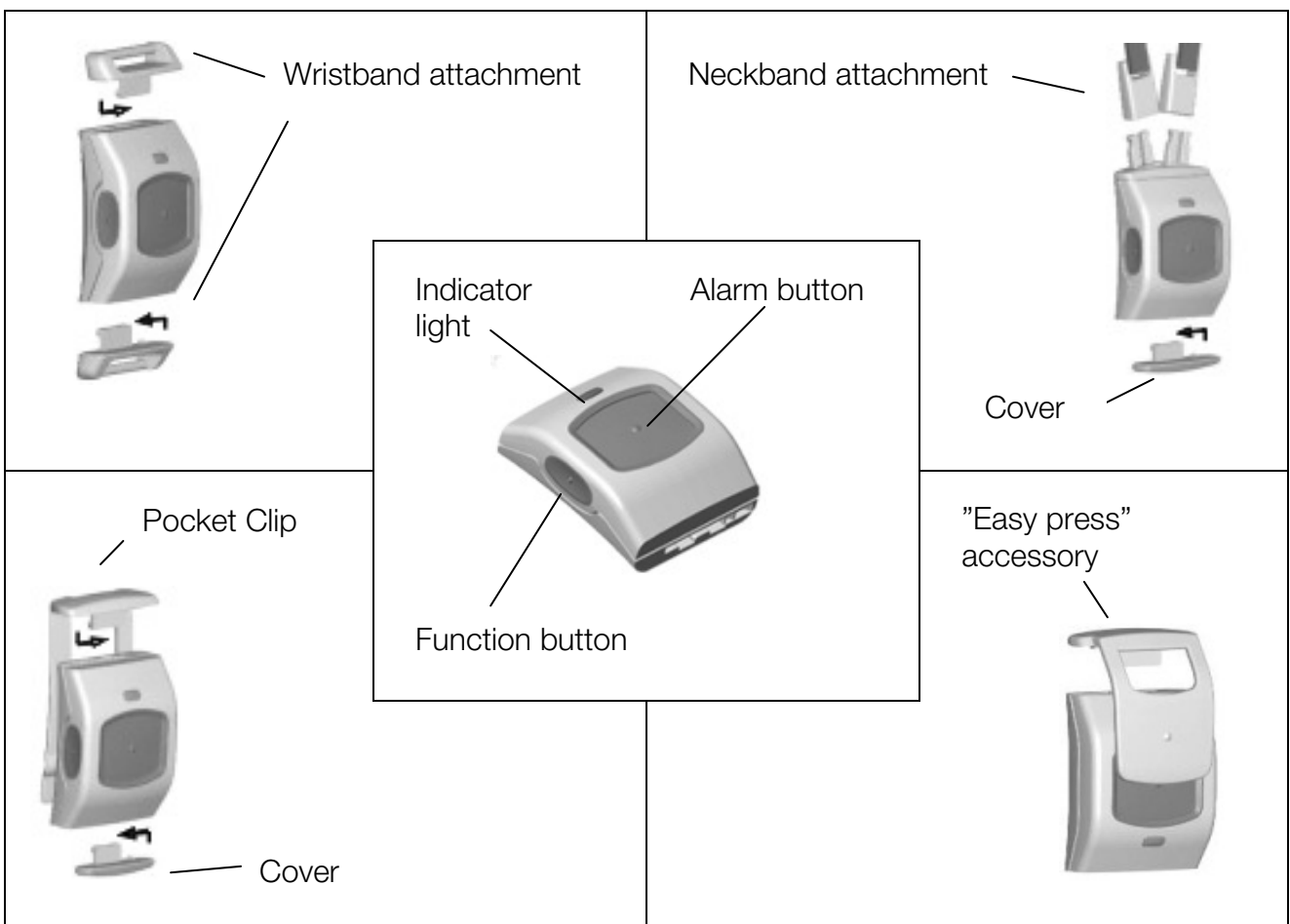
## **Description of functions**

Gina is a practical, well-designed, battery-operated carephone that is ideal for use in the individual's own home.

The Erik radio trigger is adapted for use with the Gina carephone. Erik can be worn using a wristband, neckband or pocket clip. The CareTech Easy Press Concept (EPC) accessory makes it easier for those with reduced physical strength to activate the alarm.

### **Activating the alarm**

The assistance alarm can be activated with a radio trigger or by pressing the red button on the carephone. On activation, the carephone calls an alarm receiver to summon help.

Other alarm types can be activated using accessories.

If an alarm is activated accidentally, it can be reset by pressing the green button on the carephone.

### **Different types of alarm transmitter**

Nine radio transmitters can be programmed with the carephone, e.g. a radio trigger, smoke detector, motion detector, etc.

In the event of an alarm from a radio transmitter, the battery status of the radio transmitter is tested. The Erik radio trigger has the option of automatic battery monitoring.

### **Alarm receivers**

Up to six alarm receivers can be programmed. The alarm is identified using an alarm code which is programmed into the carephone. The person receiving the alarm can then speak directly with the person raising the alarm and take appropriate action.

### **Other functions**

Additional functions can be activated or changed by programming – refer to the technical handbook or contact your supplier for further information.


## Installation

### Connection

The carephone is very simply to install. The only connection required is a connection to the telephone network via the telephone lead/adaptor.

The carephone should be connected to the master telephone socket in the home in order for it to continue to work even if a telephone receiver is left off the hook or is faulty. The master socket is normally the one closest to the entrance door – for further advice, contact your supplier.

Power is supplied via the telephone network and the built-in battery.

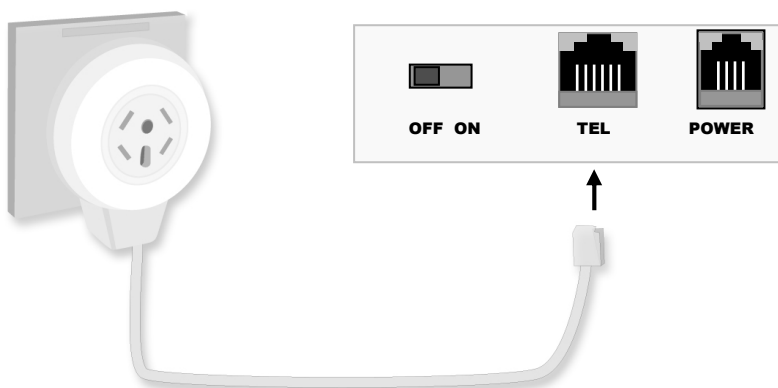- Disconnect the existing telephone from the telephone jack and connect the telephone lead/adaptor to the jack. Connect the telephone lead to the TEL socket on the carephone.
- Turn the power switch to the ON position.
- Connect the existing telephone to the intermediate plug. Lift the receiver and check that you can hear the tone.
- Carry out a test alarm from the carephone.

### Programming

The carephone's telephone numbers, alarm codes and settings are programmed via a computer, using the programming cable and the CareTech CS software.

For help with programming and information about other programming tools, please refer to the carephone's technical handbook or contact your supplier.

In some cases, the carephone can be delivered ready-programmed by the supplier.

## **Programming radio triggers**

Up to nine radio triggers or other radio-alarm transmitters, e.g. smoke detectors, can be programmed for the carephone.

Programming radio transmitters:

- Hold down the green button for approximately four seconds until the green “On” light starts to flash.
- Press the radio trigger. Once the carephone has received the radio signal, you will hear a distinct tone signal.
- Confirm by pressing the red button on the carephone.

## **Check radio coverage**

This should always be done on installation.

- Hold down the green button for approximately four seconds until the green “On” light starts to flash.
- Press the radio trigger. The built-in speaker will produce a tone signal on successful reception. Repeat this in all areas where the alarm needs to operate.

## **General information**

- If the telephone line is not working, no alarm can be sent.
- The carephone should be connected to the master telephone socket in the home.
- In order for an alarm to be sent, at least one alarm receiver/response centre needs to be programmed.
- The carephone must not come into contact with water or any other liquid.
- In case of battery alarm, the batteries should be replaced as soon as possible in order for the alarm to function correctly.
- Used batteries should be disposed of in an environmentally-responsible manner.

## **Maintenance**

### **Cleaning**

The carephone and portable alarm button should be cleaned using a slightly damp cloth. Solvents should not be used when cleaning.

### **Replacing the carephone batteries**

Batteries should be replaced within two years, regardless of whether or not the battery alarm has sounded.

- Turn off the carephone using the power switch on the reverse, disconnect the telephone lead.
- Undo the battery compartment screw.
- Open the battery compartment by removing the cover in the direction of the arrow, remove the used batteries.
- Insert new batteries, observing the correct polarity, use only LR6/AA 1,5V Alkaline batteries.
- Replace the battery compartment cover and the screws.
- Connect the telephone lead and turn the power switch to the ON position.
- Check that the unit is operating.

### **Replacing the radio trigger batteries**

To replace the radio trigger batteries, please refer to the technical handbook or contact your supplier.

## Technical data

### Carephone

Dimensions:	121 x 214 x 36 mm (LxWxH)
Power supply:	Via telephone line and built-in batteries
Batteries:	4 x 1.5 V alkaline LR6/AA
Battery life:	Up to two years
Communication:	Analogue telephone line
Radio frequency:	868.35 MHz
Radio receiver class:	3*
Number of radio transmitters:	Up to 9

### Gina

### Portable alarm button

Dimensions:	45 x 30 x 13 mm (LxWxH)
Battery:	3 V lithium battery CR2430
Battery life:	Up to five years
Water resistance:	Complies with IP67

### Erik

**Equipment class:** Class 1 radio equipment 

## Environment

This product is intended for indoor use in a normal residential environment.

Temperature:	Operating temperature +5°C to +35°C
Humidity:	0% to 75% relative humidity (non-condensed)
Environmental class:	1

\*Higher classification does not apply to this type of broadband receiver. A lower classification may possibly, but not necessarily, mean that greater attention should be paid to the risk of radio interference. Radio receiver data:

Average usable sensitivity: -100dBm, Spurious response rejection: -80dB, desensitisation: 23.3dB at 1MHz, Max usable sensitivity: N/A, Co-channel rejection: N/A

Hereby, CareTech AB, Kalix Sweden declare that this radio- and telecommunications terminal equipment is in compliance with the essential requirements and other provisions of R&TTE Directive 1999/5/EC.

The comprehensive declaration of conformity is available at the address: CareTech AB, Box 10030, SE-952 27 Kalix, Sweden

### **Environmental information**

This product complies with the requirements of the EU directive, 2002/96/EC (WEEE).

WEEE regulates the product liability for electrical and electronic recycling with the purpose of increasing recycling and minimising waste. The unit is marked with the "crossed out wheeled bin" logo, which indicates that it shall be handed in for recycling.

The product can be returned free of charge to a recycling station that is connected, directly or via a recycling system, to Care Tech or to your distributor. For detailed instructions, please check with your distributor or visit our website, [www.caretech.se](http://www.caretech.se)

Note! The WEEE information and recycling instructions applies to European Union member states only. For other countries please check local legislation or contact your distributor.

Manufactured in accordance with the EU directive, 2002/95/EC (RoHS)

The materials used in the alarm button wristband and neckband meet the textile safety requirements of Oeko-Tex standard 100.

## **Important Information**

All systems using radio and telecommunications are subject to interference beyond the user's control.

Products from CareTech are designed to minimise the impact of such interference. Nevertheless, the user must be aware that system components can be subjected to interference or other influences that may cause malfunction.

It is therefore important to regularly check that every part of the system works in all areas, especially radio communications. Contact your supplier immediately in case of any suspected malfunction.

Strong radio transmitters, such as GSM equipment, DECT telephones or wireless hearing aids, should be positioned at least two metres away from the product, since these can affect its operation. Users should pay particular attention to the risk of disruption from products which communicate using the same or adjacent frequencies.

For further information, please contact your supplier or refer to the technical handbook or the documentation at [www.caretech.se](http://www.caretech.se)

The information contained here is subject to change without notice. The only warranties for CareTech products and services are set forth in the express warranty statements accompanying such products and services. Nothing herein should be construed as constituting an additional warranty. CareTech shall not be liable for technical or editorial errors or omissions contained herein.

[www.caretech.se](http://www.caretech.se)