Grace

Bruksanvisning | User Guide

tryggitel as

Tryggitel AS
Hans Nielsen Hauges gate 50
0481 Oslo
Telefon: 22 15 75 90
firmapost@tryggitel.no
www.tryggitel.no

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
Package contents

- Grace telecare terminal
- Erik radio trigger with accessories.
- Telephone lead.
- User guide

Telecare terminal parts


There are three buttons on the front, with the following standard functions:

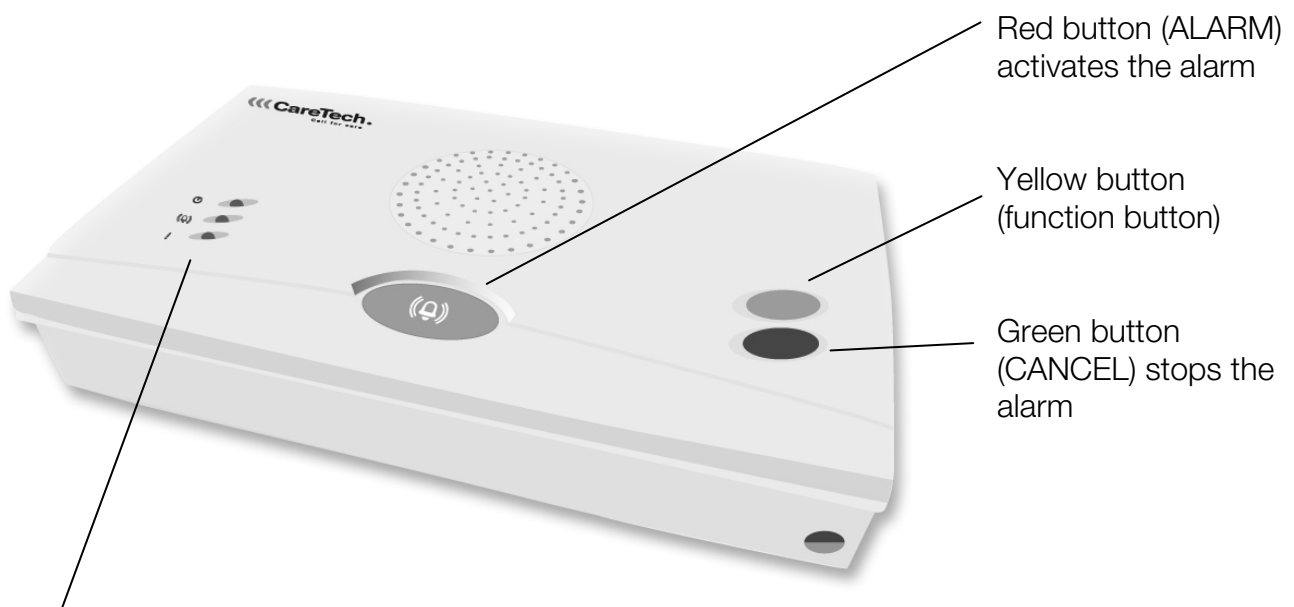The following indicator lights can be found on the front of the telecare terminal:


- ⏻ On (Green) – normally lit during operation**
 - Steady light: Indicates On state, receipt of radio signals on programming radio alarm.
 - Flashing light: Radio programming state.
- 🔔 Alarm (Red)**
 - Slow flashing: Alarm.
 - Rapid flashing: Connection to alarm receiver on alarm.
- ! Error (Red)**
 - Rapid flashing: No connection to telephone line.

The following connectors can be found on the reverse of the telecare terminal.

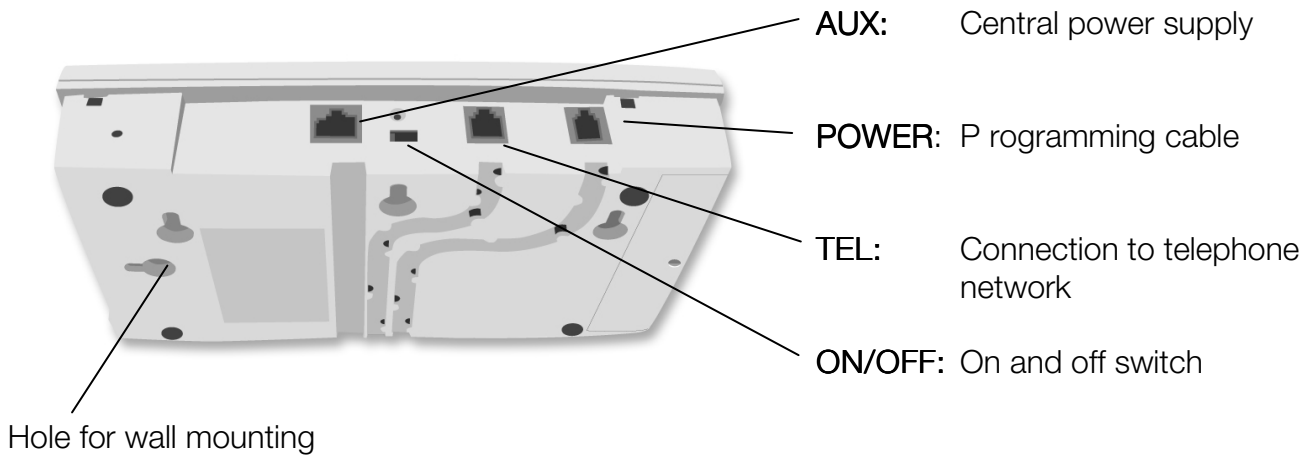Radio trigger components

Accessories included: wristband, wristband attachment, neckband, pocket clip, "easy-press" accessory and cover. The accessories are attached to the alarm button with a simple click-in action.

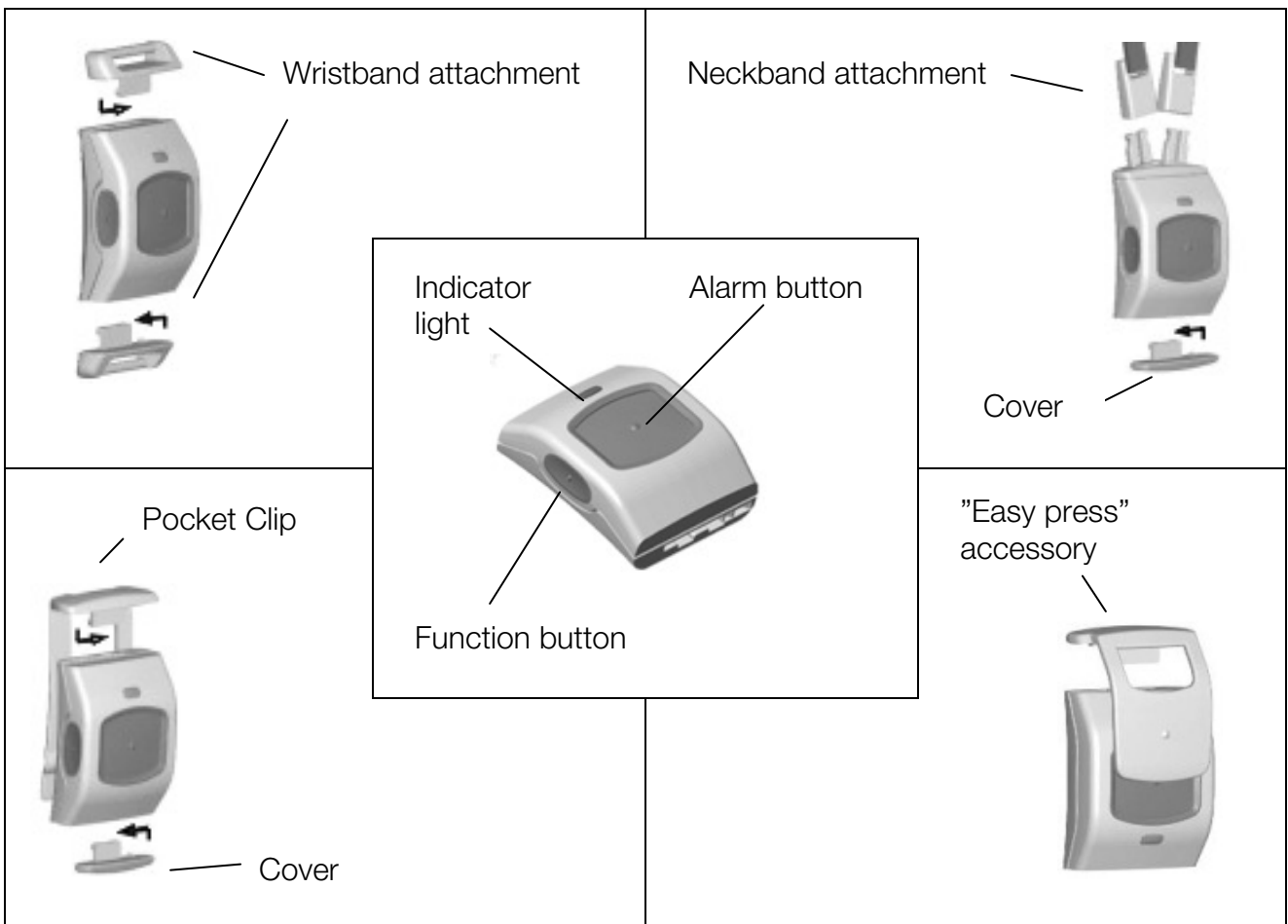Description of functions

Grace is a practical, well-designed Telecare terminal that is ideal for use in sheltered accommodation with a central power supply and telephone exchange. The Erik radio trigger is adapted for use with the Grace telecare terminal. Erik can be worn using a wristband, neckband or pocket clip. The CareTech Easy Press Concept (EPC) accessory makes it easier for those with reduced physical strength to activate the alarm.

Activating the alarm

The assistance alarm can be activated with a radio trigger or by pressing the red button on the telecare terminal. On activation, the telecare terminal calls an alarm receiver to summon help.

Other alarm types can be activated using accessories by radio, or by the wired inputs on the Telecare terminals AUX connector .

If an alarm is activated accidentally, it can be reset by pressing the green button on the telecare terminal.

Different types of alarm triggers

Nine radio transmitters can be programmed with the telecare terminal, e.g. a radio trigger, smoke detector, motion detector, etc.

Up to four wired triggers can be connected to the Telecare terminal's AUX connector, e.g. WC alarm, bed sensor, etc.

In the event of an alarm from a radio transmitter, the battery status of the radio transmitter is tested. The Erik radio trigger has the option of automatic battery monitoring.

Alarm receivers

Up to six alarm receivers can be programmed. The alarm is identified using an alarm code which is programmed into the telecare terminal. An alarm placed to a private receiver is identified by a recorded voice ID. The person receiving the alarm can then speak directly with the person raising the alarm and take appropriate action.

Other functions

Additional functions can be activated or changed by programming – refer to the technical handbook or contact your supplier for further information.


Installation

Connection

The telecare terminal is connected to the telephone network via the telephone lead/adaptor.

Central power supply and wired alarm triggers are connected to the AUX connector.

Normally the installation is made by a certified installer, for further assistance regarding installation please refer to the technical handbook or contact your supplier.

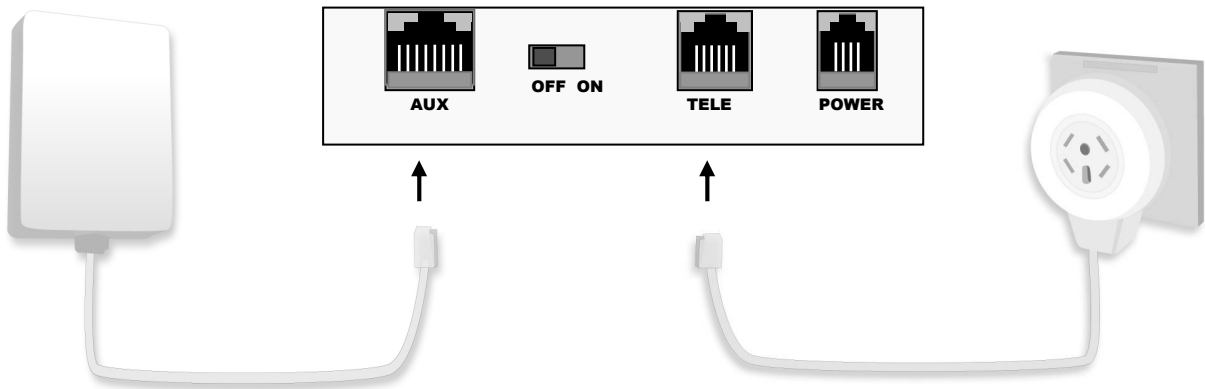Programming

The telecare terminal's telephone numbers, alarm codes and settings are programmed via a computer, using the programming cable and the CareTech CS software.

For help with programming and information about other programming tools, please refer to the telecare terminal's technical handbook or contact your supplier.

In some cases, the telecare terminal can be delivered ready-programmed by the supplier.

Programming radio triggers

Up to nine radio triggers or other radio-alarm transmitters, e.g. smoke detectors, can be programmed for the telecare terminal.

Programming radio transmitters:

- Hold down the green button for approximately four seconds until the green “On” light starts to flash.
- Press the radio trigger. Once the telecare terminal has received the radio signal, you will hear a distinct tone signal.
- Confirm by pressing the red button on the telecare terminal.

Check radio coverage

This should always be done on installation.

- Hold down the green button for approximately four seconds until the green “On” light starts to flash.
- Press the radio trigger. The built-in speaker will produce a tone signal on successful reception. Repeat this in all areas where the alarm needs to operate.

General information

- If the telephone line is not working, no alarm can be sent.
- In order for an alarm to be sent, at least one alarm receiver/response centre needs to be programmed.
- The telecare terminal must not come into contact with water or any other liquid.

Maintenance

Cleaning

The telecare terminal and portable alarm button should be cleaned using a slightly damp cloth. Solvents should not be used when cleaning.

Replacing the radio trigger batteries

To replace the radio trigger batteries, please refer to the technical handbook or contact your supplier.

Technical data

Telecare terminal

Dimensions:	121 x 214 x 36 mm (LxWxH)
Power supply:	12VDC central power supply
Inputs:	Four optoisolated
Output:	Relay output, potential-free make, max 2A
Communication:	Analogue telephone line
Radio frequency:	868.35 MHz
Radio receiver class:	3*
Number of radio transmitters:	Up to 9

Grace

Portable alarm button

Dimensions:	45 x 30 x 13 mm (LxWxH)
Battery:	3 V lithium battery CR2430
Battery life:	Up to five years
Water resistance:	Complies with IP67

Erik

Equipment class: Class 1 radio equipment 

Environment

This product is intended for indoor use in a normal residential environment.

Temperature:	Operating temperature +5°C to +35°C
Humidity:	0% to 75% relative humidity (non-condensed)
Environmental class:	1

Hereby, CareTech AB, Kalix Sweden declare that this radio- and telecommunications terminal equipment is in compliance with the essential requirements and other provisions of R&TTE Directive 1999/5/EC.

The comprehensive declaration of conformity is available at the address:
CareTech AB, Box 10030, SE-952 27 Kalix, Sweden

*Higher classification does not apply to this type of broadband receiver. A lower classification may possibly, but not necessarily, mean that greater attention should be paid to the risk of radio interference. Radio receiver data: Average usable sensitivity: -100dBm, Spurious response rejection: -80dB, desensitisation: 23.3dB at 1MHz, Max usable sensitivity: N/A, Co-channel rejection: N/A

Environmental information

This product complies with the requirements of the EU directive, 2002/96/EC (WEEE).

WEEE regulates the product liability for electrical and electronic recycling with the purpose of increasing recycling and minimising waste. The unit is marked with the "crossed out wheeled bin" logo, which indicates that it shall be handed in for recycling.

The product can be returned free of charge to a recycling station that is connected, directly or via a recycling system, to Care Tech or to your distributor. For detailed instructions, please check with your distributor or visit our website, www.caretech.se

Note! The WEEE information and recycling instructions applies to European Union member states only. For other countries please check local legislation or contact your distributor.

Manufactured in accordance with the EU directive, 2002/95/EC (RoHS)

The materials used in the alarm trigger wristband and neckband meet the textile safety requirements of Oeko-Tex standard 100.

Important Information

All systems using radio and telecommunications are subject to interference beyond the user's control.

Products from CareTech are designed to minimise the impact of such interference. Nevertheless, the user must be aware that system components can be subjected to interference or other influences that may cause malfunction.

It is therefore important to regularly check that every part of the system works in all areas, especially radio communications. Contact your supplier immediately in case of any suspected malfunction.

Strong radio transmitters, such as GSM equipment, DECT telephones or wireless hearing aids, should be positioned at least two metres away from the product, since these can affect its operation. Users should pay particular attention to the risk of disruption from products which communicate using the same or adjacent frequencies.

For further information, please contact your supplier or refer to the technical handbook or the documentation at www.caretech.se

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en

The information contained here is subject to change without notice. The only warranties for CareTech products and services are set forth in the express warranty statements accompanying such products and services. Nothing herein should be construed as constituting an additional warranty. CareTech shall not be liable for technical or editorial errors or omissions contained herein.

www.caretech.se